

# enjoy your daily class II restorations

sectional  
matrix  
systems

myClip 2.0

myQuickmat  
Forte kit

myQuickmat  
Classico kit

myCustom  
Rings kit

*powered by Style Italiano*

**INNOVATIVE TECHNIQUE**



**polydentia**  
swiss manufacture 



[www.polydentia.ch](http://www.polydentia.ch)

# myClip 2.0



All in one  
sectional matrix ring  
with integrated forceps

## Testimonial




*"myClip 2.0 is very easy to handle. Its tines have a special geometry and rigidity allowing very efficient matrix positioning on the palatine and buccal walls. This saves time for the finishing steps and ensures a good proximal morphology."*

Dr. David Gerdolle,  
Montreux - Switzerland



Literature

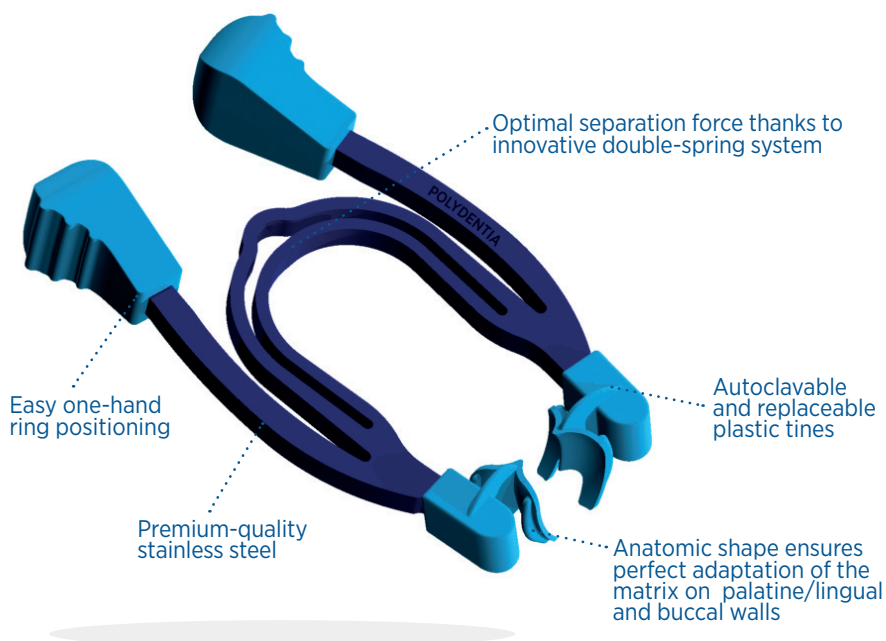

## Tines

Autoclavable and replaceable plastic tines



## LumiContrast

Dark-blue matrices which significantly increase contrast and reduce glare effect, especially while using loupes or microscopes



## Wood Wedges

Interdental wood wedges made of maple wood from sustainably managed Swiss forests



All in one sectional matrix ring with integrated forceps that allows for one-hand placement, simply using the fingers – no special ring forceps required. Its autoclavable plastic tines ensure perfect adaptation of the matrix on palatine/lingual and buccal walls, strongly reducing the need for finishing steps and therefore saving time. These anatomic shaped tines are replaceable, allowing to benefit from the long lifetime of the premium-quality stainless steel ring. The separation force of myClip 2.0 allows to create excellent contact points in the restoration of class II cavities, making it suitable for all clinical conditions.

### myClip 2.0 - REF 6305

- 1 pc. myClip 2.0,
- 2 pairs tines

### Included: Restorative trial set

- 5 pcs. LumiContrast Sectional Matrices, Premolar 0.04mm / 5mm,
- 5 pcs. LumiContrast Sectional Matrices, Molar 0.04mm / 6.4mm,
- 10 pcs. Pinky Wood Wedges, 11mm (XS),
- 10 pcs. Wood Wedges, 12mm (S)



# myQuickmat Forte kit



Complete sectional matrix system developed to offer **strong** tooth separation

## Testimonial




*"The myQuickmat Forte kit provides a strong separating ring for an effective teeth separation. The tines of myRing Forte carefully embrace the matrix and, in combination with the myWedge, ensure a perfect adaption of the matrix and a strong contact point restoration"*

Dr. Giuseppe Chiodera,  
Brescia - Italy



Literature

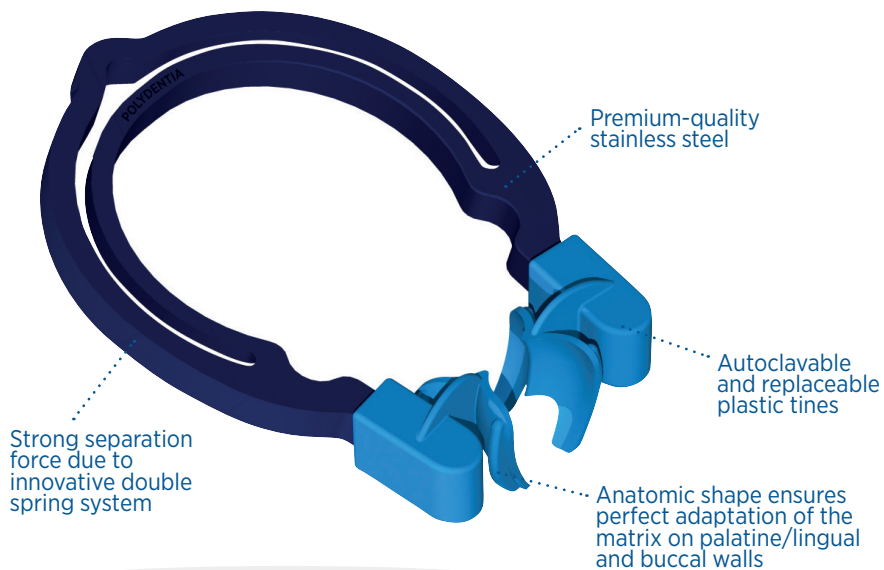

## Tines

Autoclavable and replaceable plastic tines



## LumiContrast

Dark-blue matrices which significantly increase contrast and reduce glare effect, especially while using loupes or microscopes



## myWedge

V-shaped wedge made of hi-tech plastic material. Perfectly compresses and adapts to the interproximal anatomy of teeth



## myForceps

myForceps is the indispensable instrument to safely place and remove all Polydentia sectional matrix rings



This kit features Polydentia's myRing Forte, a strong separating ring that offers optimal separation force. Its autoclavable and replaceable plastic tines ensure a perfect adaptation of the matrix to palatine/lingual and buccal walls. The kit comprises a selection of myWedge, high-quality, anatomical interdental plastic wedges. Its hollow V-shape allows it to compress and adapt to the interproximal anatomy of teeth. LumiContrast sectional matrices complete the kit's matrix assortment. These dark-blue matrices significantly increase contrast and reduce glare effect, especially while using loupes or microscopes.

### myQuickmat Forte kit - REF 6803

- 50 pcs. Quickmat Sectional Matrices, Premolar 0.04mm / 5mm,
- 50 pcs. Quickmat Sectional Matrices, Molar 0.04mm / 6.4mm,
- 20 pcs. Quickmat Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.,
- 20 pcs. LumiContrast Sectional Matrices, Premolar 0.04mm / 5mm,
- 20 pcs. LumiContrast Sectional Matrices, Molar 0.04mm / 6.4mm,
- 20 pcs. LumiContrast Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.,
- 2 pcs. myRing Forte,
- 2 pairs tines,
- 50 pcs. myWedge, Extrasmall (XS),
- 50 pcs. myWedge, Small (S),
- 40 pcs. myWedge, Medium (M),
- 40 pcs. myWedge, Large (L),
- 1 pc. myForceps



# myQuickmat Classico kit




Complete sectional matrix system developed to offer **gentle** tooth separation

## Testimonial



*"myRing Classico applies the perfect amount of pressure to ensure optimal matrix adaptation and retention. It does not lose its shape: its mechanical performance is consistent and does not diminish with use. The newly re-designed Delta Tubes - powered by Style Italiano - provide significant versatility to perfectly adapt the sectional matrix according to the case at hand."*


Dr. Jordi Manauta,  
Italy

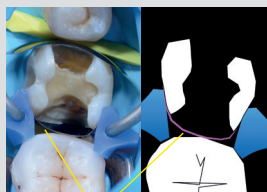

"No matrix deformation even in case of broad cavities"



Literature

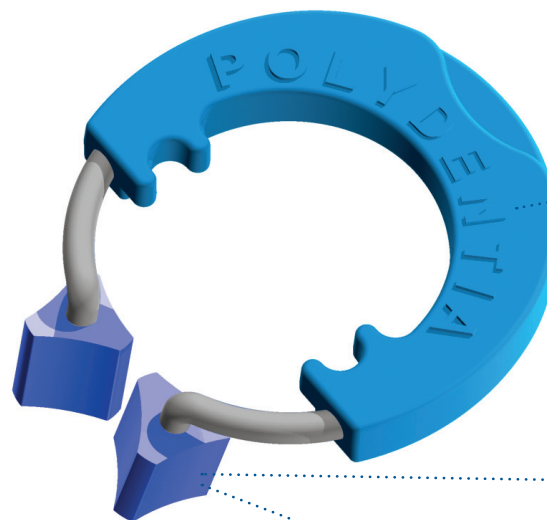

Nickel-Titanium ring with reinforcement

Autoclavable and replaceable Delta Tubes

Triangular silicone tubes allow perfect adaptation of the matrix

## Delta Tubes

Anatomically shaped silicone tubes which ensure optimal matrix retention and adaptation to the tooth's anatomy



## Steel Matrices

For cervical, subgingival restorations



## LumiContrast

Dark-blue matrices which significantly increase contrast and reduce glare effect, especially while using loupes or microscopes



## Wood Wedges

Interdental wood wedges made of maple wood from sustainably managed Swiss forests



This kit features Polydentia's myRing Classico, a classic sectional matrix ring made of nickel-titanium. It is reinforced to provide an optimal and gentle separation force in combination with the included Interdental Maple Wood Wedges. The kit also includes the newly re-designed Delta Tubes. These anatomically shaped silicone tubes ensure optimal matrix retention and adaptation to the tooth's anatomy. LumiContrast sectional matrices complete the kit's assortment. These dark-blue matrices significantly increase contrast and reduce glare effect, especially while using loupes or microscopes.

### myQuickmat Classico kit - REF 6801

- 50 pcs. Quickmat Sectional Matrices, Premolar 0.04mm / 5mm,
- 50 pcs. Quickmat Sectional Matrices, Molar 0.04mm / 6.4mm,
- 20 pcs. Quickmat Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.,
- 20 pcs. LumiContrast Sectional Matrices, Premolar 0.04mm / 5mm,
- 20 pcs. LumiContrast Sectional Matrices, Molar 0.04mm / 6.4mm,
- 20 pcs. LumiContrast Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.,
- 2 pcs. myRing Classico,
- 6 pairs Delta Tubes,
- 50 pcs. Pinky Wood Wedges, 11mm (XS),
- 50 pcs. Wood Wedges, 12mm (S),
- 40 pcs. Wood Wedges, 13mm (M),
- 40 pcs. Wood Wedges, 15mm (L)



# myCustom Rings Kit

INNOVATIVE TECHNIQUE

powered by  
STYLE ITALIANO




## Testimonial




*"In almost any case where the proximal embrasures are healthy and intact, and there is the risk of losing them, myCustom Rings is the best technique to overcome this common difficulty. Easy to use, as well as easy to produce, this is a strategic technique for the everyday adhesive modern dentistry."*

Dr. Anna Salat,  
Italy

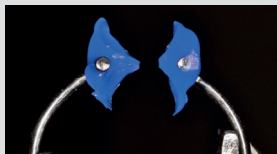


or  
[www.styleitaliano.org/  
mycustom-rings-kit](http://www.styleitaliano.org/mycustom-rings-kit)

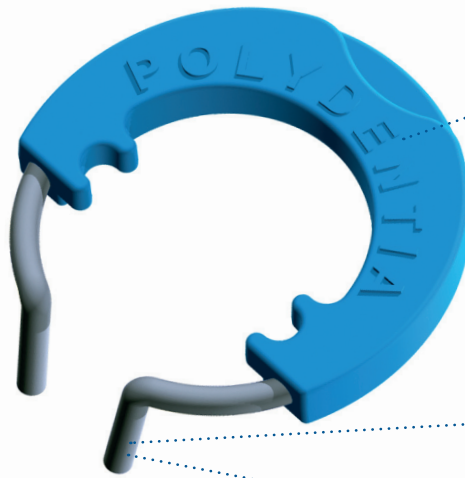

Nickel-Titanium ring with reinforcement

Custom shaped tines made using light curable myCustom Resin

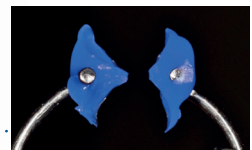

Custom tines enable to achieve strong contact points and precise proximal contour during restoration

## myCustom Resin

The ring tines are custom-shaped by taking an impression of the healthy and intact proximal walls using the included light curable myCustom Resin




## LumiContrast

Dark-blue matrices which significantly increase contrast and reduce glare effect, especially while using loupes or microscopes



## Wood Wedges

Interdental wood wedges made of maple wood from sustainably managed Swiss forests

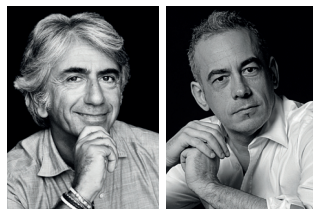

"The innovative technique behind myCustom Rings kit perfectly exemplifies our philosophy 'Feasible, teachable and repeatable: the Style Italiano recipe for a daily success dentistry'."

**Prof. Angelo Putignano MD, DDS  
& Dr. Walter Devoto, DDS  
Founders of Style Italiano -Italy**

myCustom Rings kit has been developed in partnership with Style Italiano, a study group of clinicians and researchers specialized in aesthetic and restorative dentistry. This kit is conceived for the application of an innovative technique developed by Dr. Jordi Manauta of Style Italiano. The myCustom Rings technique allows to easily recreate the pre-existing interproximal anatomy. The ring tines are custom-shaped by taking an impression of the healthy and intact proximal walls using the included light curable myCustom Resin. This enables to achieve strong contact points and precise proximal contour during restoration.

### myCustom Rings kit - REF 6802


- 50 pcs. Quickmat Sectional Matrices, Premolar 0.025mm / 5mm,
- 50 pcs. Quickmat Sectional Matrices, Molar 0.025mm / 6.4mm,
- 20 pcs. LumiContrast Sectional Matrices, Premolar 0.025mm / 5mm,
- 20 pcs. LumiContrast Sectional Matrices, Molar 0.025mm / 6.4mm,
- 20 pcs. Quickmat Transparent Sectional Matrices, Premolar 0.075mm / 5mm,
- 20 pcs. Quickmat Transparent Sectional Matrices, Molar 0.075mm / 6.4mm,
- 2 pcs. myRing Classico,
- 6 pairs Delta Tubes,
- 50 pcs. Pinky Wood Wedges, 11mm (XS), 50 pcs. Wood Wedges, 12mm (S),
- 40 pcs. Wood Wedges, 13mm (M), 40 pcs. Wood Wedges, 15mm (L),
- 1 pc. myCustom Resin, 3ml syringe with 10 pcs. needle tips gauge 18 for myCustom Resin



# myWedge

myWedge is an interproximal wedge made of hi-tech plastic material. Its innovative hollow v-shape allows it to compress and adapt to the interproximal anatomy of teeth. It prevents damage to the papilla while guiding and forming the matrix band to the anatomy of the cervical area of the tooth. myWedge's unique round-head geometry has been designed for easy handling, placement, and removal of the wedge with any standard tweezers or forceps.

- **Innovative hollow v-shape**
- **Gripped head to prevent sliding**
- **Flashy colors for easy identification**

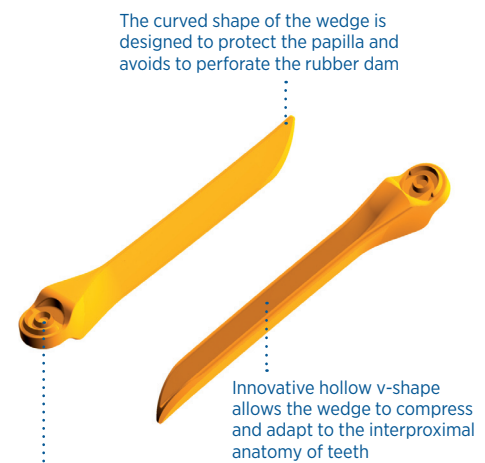


myWedge's unique round-head geometry has been designed for easy handling, placement, and removal of the wedge with any standard tweezers or forceps

# LumiContrast

LumiContrast sectional matrices are dark-blue matrices which significantly increase contrast and reduce glare effect, especially while using loupes or microscopes.

- **High contrast**
- **Low reflection**
- **Ideal when using loupes or microscopes**

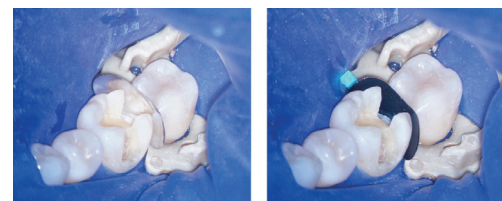

Without LumiContrast

With LumiContrast

# myQuickmatrix Forceps

myQuickmatrix Forceps ensures a firm grip of the matrix when placing and removing it during restoration.

- **For easy matrix placement and removal**
- **Tips with special ribbing for a firm grip**
- **Resistant tips in tungsten carbide**



# myForceps

myForceps is the indispensable instrument to safely place and remove all Polydentia sectional matrix rings such as myRing Classico and myRing Forte. The shaped tips ensure a firm grip of the ring during application.

- **For easy placement of myRing Classico and myRing Forte**
- **Strong and resistant**
- **Outstanding quality**



# Kit refills

## Sectional Matrix Rings



**REF 6805**  
myClip 2.0  
2 pcs.,  
Tines  
2 pairs



**REF 6804**  
myRing  
Classico  
2 pcs.,  
Delta Tubes  
2 pairs



**REF 6806**  
myRing  
Forte  
2 pcs.,  
Tines  
2 pairs









**REF 6303**  
8 pairs  
Tines for  
myClip 2.0  
& myRing  
Forte






**REF 5712**  
15 pairs  
Delta Tubes









## Steel Matrices

	Quantity	REF
 Quickmat Sectional Matrices, Premolar 0.025mm / 5mm	100	5735
 Quickmat Sectional Matrices, Molar 0.025mm / 6.4mm	100	5736
 Quickmat Sectional Matrices, Premolar 0.04mm / 5mm	100	5705
 Quickmat Sectional Matrices, Molar 0.04mm / 6.4mm	100	5706
 Quickmat Sectional Matrices, Extended 0.025mm / 6.4mm+Ext.	50	5756
 Quickmat Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.	50	5746









## Transparent Matrices

	Quantity	REF
 Quickmat Transparent Sectional Matrices, Premolar 0.075mm / 5mm	100	5765
 Quickmat Transparent Sectional Matrices, Molar 0.075mm / 6.4mm	100	5766
 Quickmat Transparent Sectional Matrices, Extended 0.075mm / 6.4mm+Ext.	50	5776

## LumiContrast Matrices

	Quantity	REF
 LumiContrast Sectional Matrices, Premolar 0.025mm / 5mm	70	6735
 LumiContrast Sectional Matrices, Molar 0.025mm / 6.4mm	70	6736
 LumiContrast Sectional Matrices, Premolar 0.04mm / 5mm	70	6705
 LumiContrast Sectional Matrices, Molar 0.04mm / 6.4mm	70	6706
 LumiContrast Sectional Matrices, Extended 0.025mm / 6.4mm+Ext.	50	6756
 LumiContrast Sectional Matrices, Extended 0.04mm / 6.4mm+Ext.	50	6746

## Wedges

	Quantity	REF
 Pinky Wood Wedges, 11mm (XS)	200	5011
 Wood Wedges, 12mm (S)	100	5012
 Wood Wedges, 13mm (M)	100	5013
 Wood Wedges, 15mm (L)	100	5015
 myWedge, Extrasmall (XS)	100	7101
 myWedge, Small (S)	100	7102
 myWedge, Medium (M)	100	7103
 myWedge, Large (L)	100	7104

## Resin

	REF
myCustom Resin, 3ml Syringe & 10 pcs. needle tips gauge 18	6811

## Forceps

	REF
myQuickmatrix Forceps	5703
myForceps	5707

# let's keep connected



REV 05 03.2018



Polydentia SA  
via Cantonale 47  
6805 Mezzovico-Vira  
Switzerland

phone: +41 91 946 29 48  
fax: +41 91 946 32 03

info@polydentia.ch  
www.polydentia.ch

